

HIP MODULARITY

REMEDY® MODULAR STEMS


The REMEDY® Hip Spacer is part of the treatment foreseen in a two-stage procedure performed in the event of permanent prosthesis infection. The REMEDY® Hip Spacer implant is intended for temporary use only (180 days or less). It allows basic joint mobility and releases antibiotics into the joint area to protect the implant from bacterial colonization. A second surgery will be required at a later date to remove the REMEDY® Hip Spacer and replace it with a permanent hip joint implant.



Stainless steel rod reinforces the femoral stem

REMEDY® Hip Spacers:

- Single-use medical devices/ethylene oxide sterile
- Formed with bone cement (PMMA) and gentamicin

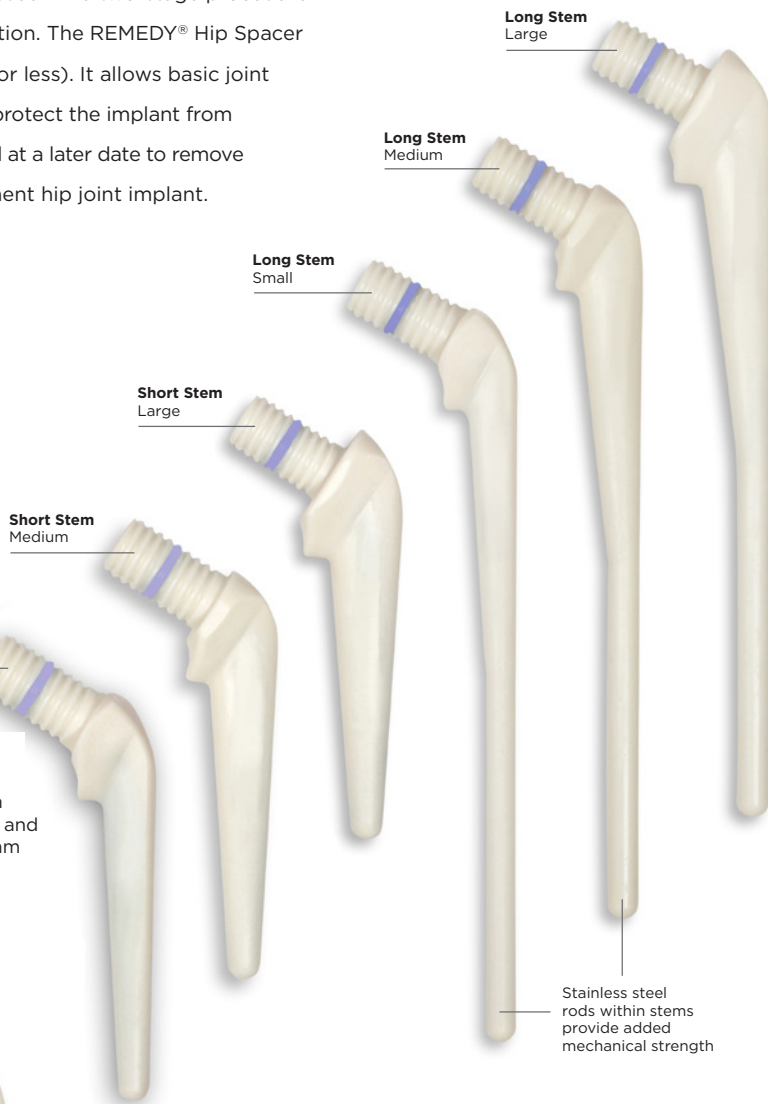

REMEDY® ACETABULAR CUP

- The first available pre-formed PMMA, antibiotic eluting acetabular cup spacer
- Cup articulates with head to improve device stability and mitigate dislocations
- Helps protect acetabular integrity during staged procedures

REMEDY® MODULAR HEADS



*40mm head size only available in REMEDY SPECTRUM® GV System

OsteoRemedies, LLC | 6800 Poplar Avenue | Suite 120 | Memphis, TN 38138
 1-800-OSTEO-XL | 901-453-3141 | info@OsteoRemedies.com | OsteoRemedies.com

