

R E S O R B A B L E   B E A D   K I T   ( R B K )

# OSTEOBOOST<sup>®</sup> SELECT<sup>™</sup>

100% PURE, SYNTHETIC,  
RESORBABLE BONE VOID FILLER  
WITH NANO-CRYSTALLINE TECHNOLOGY

OVER 20,000 CLINICAL CASES



Highly Moldable  
Putty Packs Into Voids  
and can be Used to Coat  
Implant Surfaces

**Bioengineered Triphasic Graft,  
Enhanced with nHA**


- Calcium Sulfate
- Beta-tricalcium Phosphate
- Nano-crystalline Hydroxyapatite (nHA)

Biocompatible


Neutral pH

Optimized Resorption Rate

Indicated for  
Gaps and Voids Located  
in Extremities, Spine,  
Pelvis, or Cranium

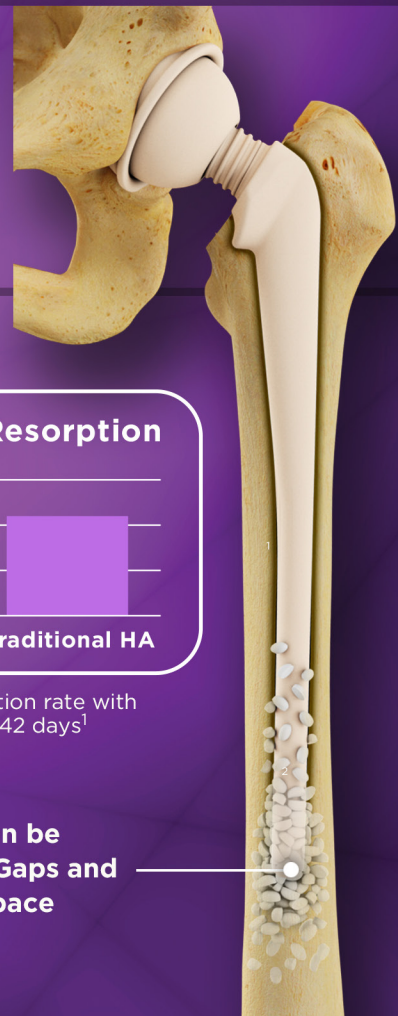

## Months to Resorption

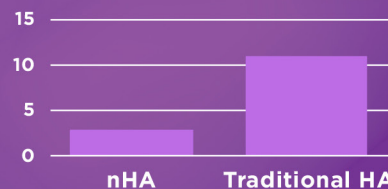

Optimized resorption rate with  
dissolution in 42 days<sup>1</sup>

Beads can be  
Packed Into Gaps and  
Dead Space

# OSTEOBOOST® SELECT™ RESORBABLE BEAD KIT

*Complex Problem - Simple Solution*

**OSTEOBOOST® SELECT™ RBK** utilizes a time-tested bone void filler with over 10 years clinical use to fill voids and gaps not intrinsic to the stability of the bone structure.

## Sterile Resorbable Bead Convenience Kit

- Bone Void Filler Powder
- Distilled Water
- Mixing Spatula
- Spreading Spatula
- Bead Mat



## Elongated Bead Shape Maximizes Surface Area



**X-Small**  
(3 x 4.5mm)



**Small**  
(4.5 x 7mm)




**Medium**  
(5.5 x 9.5mm)

**OSTEOBOOST® SELECT™ RBK's** unique elongated bead shape is available in three bead sizes and one double-sided bead mat for both 5cc and 10cc kits.

**Front of Mat with Small and  
Medium Beads**

**Back of Mat with  
X-Small Beads**

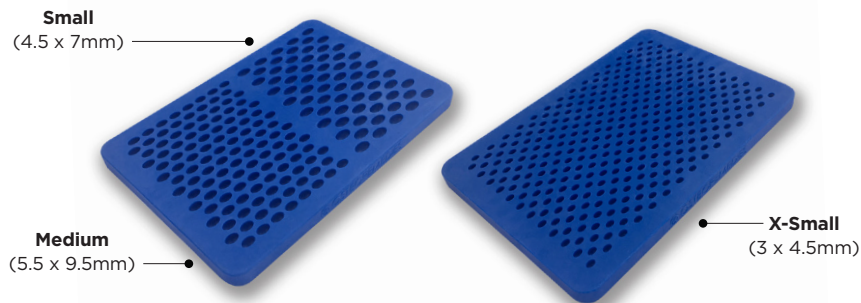

OsteoRemedies, LLC | 6800 Poplar Avenue | Suite 120 | Memphis, TN 38138

1-800-OSTEO-XL | 901-453-3141 | [info@OsteoRemedies.com](mailto:info@OsteoRemedies.com) | [OsteoRemedies.com](http://OsteoRemedies.com)

OsteoRemedies, OsteoBoost and the corporate mark are registered trademarks of OsteoRemedies, LLC  
1. Data on file. ©2024 • MKT032 Rev 05 02-2024



**OSTEOREMEDIES®**  
ADVANCED MEDICAL TECHNOLOGIES